

# Superbrush



## How to brush


Hug.



**1.** Place Superbrush over the teeth and brush a few teeth at a time

Superbrush cleans 3 sides of teeth at the same time and makes brushing more effective.

Tilt.

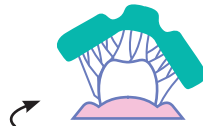

**2.** Tilt Superbrush slightly for optimal plaque removal along the gum line

Superbrush hugs the teeth to improve and simplify brushing for children and adults.

Brush!



**3.** Use a light grip and brush with short, gentle scrubbing strokes

Give your teeth a hug with superbrush!

The bristles are precisely angled in 45° to effectively clean along the gumline and prevent build up of harmful plaque.